beurer

FC 76

DE Mikrodermabrasionsgerät Gebrauchsanweisung	2	TR Mikrodermabrazyon cihazı Kullanım kılavuzu	57
EN Microdermabrasion device Instructions for use	13	RU Прибор по уходу за кожей лица Инструкция по применению.....	68
FR Appareil de microdermabrasion Mode d'emploi.....	24	PL Urządzenie do mikrodermabrazji Instrukcja obsługi.....	79
ES Aparato de microdermoabrasión Instrucciones de uso.....	35		
IT Apparecchio per microdermoabrasione Istruzioni per l'uso.....	46		

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Read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.

WARNING

- The device is only intended for domestic/private use, not for commercial use.
- This device is not suitable for persons under the age of 18.
- This device may be used by people with reduced physical, sensory or mental skills or a lack of experience or knowledge, provided that they are supervised or have been instructed on how to use the device safely and are fully aware of the consequent risks of use.
- Children must not play with the device.
- Cleaning and user maintenance must not be performed by children unless supervised.
- If the mains connection cable of this device is damaged, it must be disposed of. If it cannot be removed, the device must be disposed of.
- Under no circumstances should you open or repair the device yourself, as faultless functionality can no longer be guaranteed thereafter. Failure to comply will invalidate the warranty.

Dear customer,

Thank you for choosing a product from our range. Our name stands for high-quality, thoroughly tested products for applications in the areas of heat, weight, blood pressure, body temperature, pulse, gentle therapy, massage, beauty, air and baby.

With kind regards

Your Beurer team

Getting to know your instrument

The FC 76 enables professional facial peels with microdermabrasion technology. A high-quality sapphire coating is used to gently remove excess skin cells and the state-of-the-art vacuum technology noticeably improves skin texture and enhances circulation. The treatment also simulates cell renewal for a finer and smoother complexion.

Table of contents

1. Included in delivery.....	14	8. Cleaning and maintenance.....	20
2. Signs and symbols	15	9. Accessories and replacement parts	21
3. Intended use.....	16	10. What if there are problems?	21
4. Warnings and safety notes	16	11. Disposal.....	21
5. Device description.....	17	12. Technical specifications	22
6. Initial use	17	13. Warranty / Service.....	22
7. Usage	18		

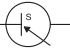
1. Included in delivery

Check that the exterior of the cardboard delivery packaging is intact and make sure that all contents are present. Before use, ensure that there is no visible damage to the device or accessories and that all packaging material has been removed. If you have any doubts, do not use the device and contact your retailer or the specified Customer Service address.

- 1x FC 76
- 1x precision attachment
- 1x coarse attachment
- 1x fine attachment
- 20x single-use filter
- 1x mains adapter
- These instructions for use

2. Signs and symbols

The following symbols are used on the device, in these instructions for use, on the packaging and on the type plate for the device:

	Warning Warning instruction indicating a risk of injury or damage to health		Manufacturer
	Important Safety note indicating possible damage to the device/ accessory		Product information Note on important information
	Danger The device must not be used near to water or in water (e.g. in a wash basin, or in the shower or bath) – risk of electric shock!		The products demonstrably meet the requirements of the Technical Regulations of the EAEU.
	Read the instructions		Dispose of packaging in an environmentally friendly manner
	Protection class II device The device is double-protected and thus corresponds to protection class 2		CE labelling This product satisfies the requirements of the applicable European and national directives.
	Disposal in accordance with the Waste Electrical and Electronic Equipment EC Directive – WEEE		Only use indoors.
IP20	Protected against solid foreign objects 12.5 mm in diameter and larger		Polarity of d.c. power connector
	Power supply unit, switch mode; switch mode power supply unit; SMPS		Safety isolating transformer, short-circuit proo

3. Intended use

The device is only to be used on the face (apart from the eye area and mouth) of adults. Do not use the device on animals!

The device is only intended for the purpose described in these instructions for use. The manufacturer is not liable for damage resulting from improper or careless use.

4. Warnings and safety notes

WARNING

- Danger of suffocation by packaging material. Keep packaging material away from children.
- Before use, inspect the device and accessories for visible damage. If you have any doubts, do not use the device and contact your retailer or the specified Customer Services address.
- Keep the mains adapter away from water.

DANGER

This device must not be used near to water or in water (e.g. in a wash basin, or in the shower or bath) – risk of electric shock!

WARNING

Do **NOT** use the device

- on cracked skin
- on open wounds
- in the eye area or on the mouth
- on any part of the body other than the face (except for eye area and mouth)
- if you have a skin disease (e.g. herpes, vascular diseases or skin growths)
- if you have skin cancer
- if you are suffering from skin irritation
- immediately after sunbathing or if you have sunburn
- if you are taking medication that contains steroids
- if you are taking blood thinners
- if you have multiple sclerosis, diabetes or an autoimmune disorder

WARNING

For this reason, the device must only be operated as follows only with the supplied power supply unit at the mains voltage indicated on the power supply unit.

IMPORTANT

Only specialist personnel may perform repairs on electrical devices. Improper repairs may subject users to considerable danger. In the event of a fault or damage, have the device repaired by a certified specialist workshop.

WARNING

- Only use the device in dry rooms (e.g. never in the bath, sauna).
- Never reach for a device that has been submerged in water. Disconnect the mains adapter immediately.
- Do not drop the device or subject it to any impacts.
- Do not clamp or pinch the mains cable.
- Ensure that the mains cable and mains adapter do not come into contact with water, steam or other liquids.

5. Device description

1. Precision attachment with fine sapphire coating, suitable for precise treatment of small areas (e.g. in nose area)
2. Fine attachment with fine sapphire coating, suitable for sensitive skin
3. Coarse attachment with coarse sapphire coating, suitable for intensive treatment of problem areas
4. Hand-held device
5. ON/OFF switch and adjustment of intensity levels 1-2

6. Initial use

Note

Before you start using the device, make sure that your face is clean, dry and free of any traces of products such as powder, cream or make-up.

6.1 Carrying out a skin test

Carry out a skin test 24 hours before using the device. Test the device at the lowest intensity level on small, inconspicuous areas of your neck. If after 24 hours the test-area is red or painful, do not use the device.

WARNING

- Do not use the device more than once a week.
- Immediately after use, the treated skin may be slightly reddened and in very rare cases there may be small spots of bleeding under the skin.
- For one day before and one day after use, avoid activities that cause you to perspire, such as saunas, tanning beds and sun exposure.


Note


The device can also be used twice a week for non-sensitive problem skin (with a break of at least 48 hours between use).

7. Usage


1. Connect the mains adapter to the device and plug the adapter into the socket.

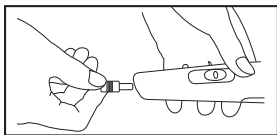
2. Insert a single-use filter into the head of the device.

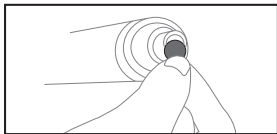3. Place the attachment you wish to use on the tip of the hand-held device.

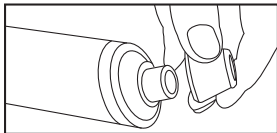4. To switch on the device, move the ON/OFF switch to the desired intensity level.

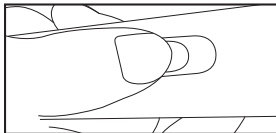
5. Place the hand-held device on your face and move it over the skin. Make sure that you do not move the hand-held device over one part of the face more than once. Further details on treatment can be found in Section "6.3 Application areas".

i Note

The pressure you apply when using the device should be comfortable and not painful. Some slight reddening should be expected after use.

6. When you have completed the application, switch off the device using the ON/OFF switch and disconnect the mains adapter from the socket..

7. Change the single-use filter after each use. To do this, pull the attachment off from the top. Use tweezers to remove the single-use filter.

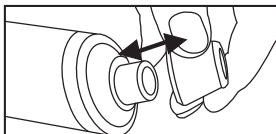
Application areas

i Note

Ensure that the part of your face you are treating is taut when you move the hand-held device over the skin. You can do this by pulling the skin taut to the side or downwards with your free hand.

Forehead:

Move the hand-held device slowly and evenly outwards from the centre of the forehead. Make sure that you do not treat any part of your skin more than once.

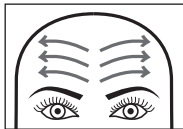Temples:

Move the hand-held device outwards towards your hairline. Repeat this to the left and right of the nose. Ensure that the device is at least 1 cm away from your eyes.

**Nose:**

Move the hand-held device downwards from the centre at the top of the bridge of your nose. Repeat this to the left and right of the nose.

**Cheeks:**

Move the hand-held device outwards across your cheeks.

**Around the mouth:**

Move the hand-held device along your upper lip from the centre to the left and right. Move the hand-held device upwards from your bottom lip to the left and right. Make sure that you do not touch your lips.

**8. Cleaning and maintenance**

- Never submerge the device in water or other liquids.
- Disinfect the attachments with a standard disinfectant after each use.
- Change the single-use filter after each use.
- Use a dry cloth to clean the device of any dirt. Do not use any abrasive cleaning products.
- Protect the device from knocks, damp, dust, chemicals, marked temperature fluctuations, and nearby sources of heat (ovens, heaters).

9. Accessories and replacement parts

For replacement parts and additional items visit www.beurer.com or contact the corresponding service address in your country (according to the service address list).

Replacement parts and additional items are also available from retailers.

Replacement part	Item number and/or order number
10 single-use filters	163.660
Precision attachment	163.767
Coarse attachment	163.766
Large attachment	163.765

10. What if there are problems?

Problem	Cause	Solution
Device will not switch on	On/Off switch faulty	Contact Customer Services or your retailer
No power	The mains part is not connected correctly.	Connect the mains part correctly.
No power	The mains part is faulty.	Contact Customer Services or your retailer
The device's suction is incorrect	Filter blocked/pump faulty	Replace the old filter with a new filter
Exfoliation effect from attachments is no longer ensured	Sapphire coating has worn down	Contact Customer Services or your retailer to request new attachments.

11. Disposal

For environmental reasons, do not dispose of the device in the household waste at the end of its useful life. Dispose of the unit at a suitable local collection or recycling point.

Dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment). If you have any questions, please contact the local authorities responsible for waste disposal.

12. Technical specifications

Power unit voltage	Input: 100–240V AC / 50/60 Hz Output: 5,0 V DC / 1,0 A / 5,0 W
Dimensions	190 x 40 x 39 mm
Weight	~128 g
Average active efficiency	≥ 77,32%
No-load power consumption	≤ 0,049 W

13. Warranty / Service

Beurer GmbH, Söflinger Straße 218, 89077 Ulm, Germany (hereinafter referred to as “Beurer”) provides a warranty for this product, subject to the requirements below and to the extent described as follows.

The warranty conditions below shall not affect the seller’s statutory warranty obligations which ensue from the sales agreement with the buyer. The warranty shall apply without prejudice to any mandatory statutory provisions on liability.

Beurer guarantees the perfect functionality and completeness of this product.

The worldwide warranty period is 3 years, commencing from the purchase of the new, unused product from the seller.

The warranty only applies to products purchased by the buyer as a consumer and used exclusively for personal purposes in the context of domestic use. German law shall apply.

During the warranty period, should this product prove to be incomplete or defective in functionality in accordance with the following provisions, Beurer shall carry out a repair or a replacement delivery free of charge, in accordance with these warranty conditions.

If the buyer wishes to make a warranty claim, they should approach their local retailer in the first instance: see the attached “International Service” list of service addresses.

The buyer will then receive further information about the processing of the warranty claim, e.g. where they can send the product and what documentation is required.

A warranty claim shall only be considered if the buyer can provide Beurer, or an authorised Beurer partner, with

- a copy of the invoice/purchase receipt, and
- the original product.

The following are explicitly excluded from this warranty:

- deterioration due to normal use or consumption of the product;
- accessories supplied with this product which are worn out or used up through proper use (e.g. batteries, rechargeable batteries, cuffs, seals, electrodes, light sources, attachments and nebuliser accessories);
- products that are used, cleaned, stored or maintained improperly and/or contrary to the provisions of the instructions for use, as well as products that have been opened, repaired or modified by the buyer or by a service centre not authorised by Beurer;
- damage that arises during transport between manufacturer and customer, or between service centre and customer;
- products purchased as seconds or as used goods;
- consequential damage arising from a fault in this product (however, in this case, claims may exist arising from product liability or other compulsory statutory liability provisions).

Repairs or an exchange in full do not extend the warranty period under any circumstances.

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